

WEGO BLOOD BAG

REDUCE PAIN, CONVENIENT BLOOD STORAGE



BLOOD COLLECTION



BLOOD TRANSFUSION



BLOOD STORAGE



SEPARATE BLOOD COMPONENTS

MULTIPLE CHOICES

Meet your needs



CALENDERING



BLOOD BAG CPDA-1 SOLUTION

STERILE, PYROGEN-FREE.
For collection of 450 ml of Blood.
63ml of Anticoagulant Citrate Phosphate Dextrose
Solution (USP) Contains:
Each 100ml of CPDA-1 Contains:
Adrenaline monohydrate 2.63g
Citric Acid dihydrate 5.225g
Sodium Citrate dihydrate 3.15g
Sodium Biphosphate monohydrate 0.2275g
Dextrose monohydrate @ 5

CAUTION

1. Do not use unless the anticoagulant is clear.
2. Store blood between + 2°C and + 6°C.
3. Crossmatch before transfusion.
4. Mix blood thoroughly immediately before use.
5. Do not add medication to this blood.
6. Transfusion set must have a filter.

DO NOT VENT

Batch No: 20220410
Mfg. Date: 20240410
Exp. Date:

Donor Number (Name)

Blood Group Rh - Type

Collection Date

Do not infuse after

Serology: non-reactive to

Test




BLOW MOLDING

DUAL PACKAGING

Eliminates microbial contamination on the surface
and maintains the content of anticoagulant

ALUMINUM PACKAGING

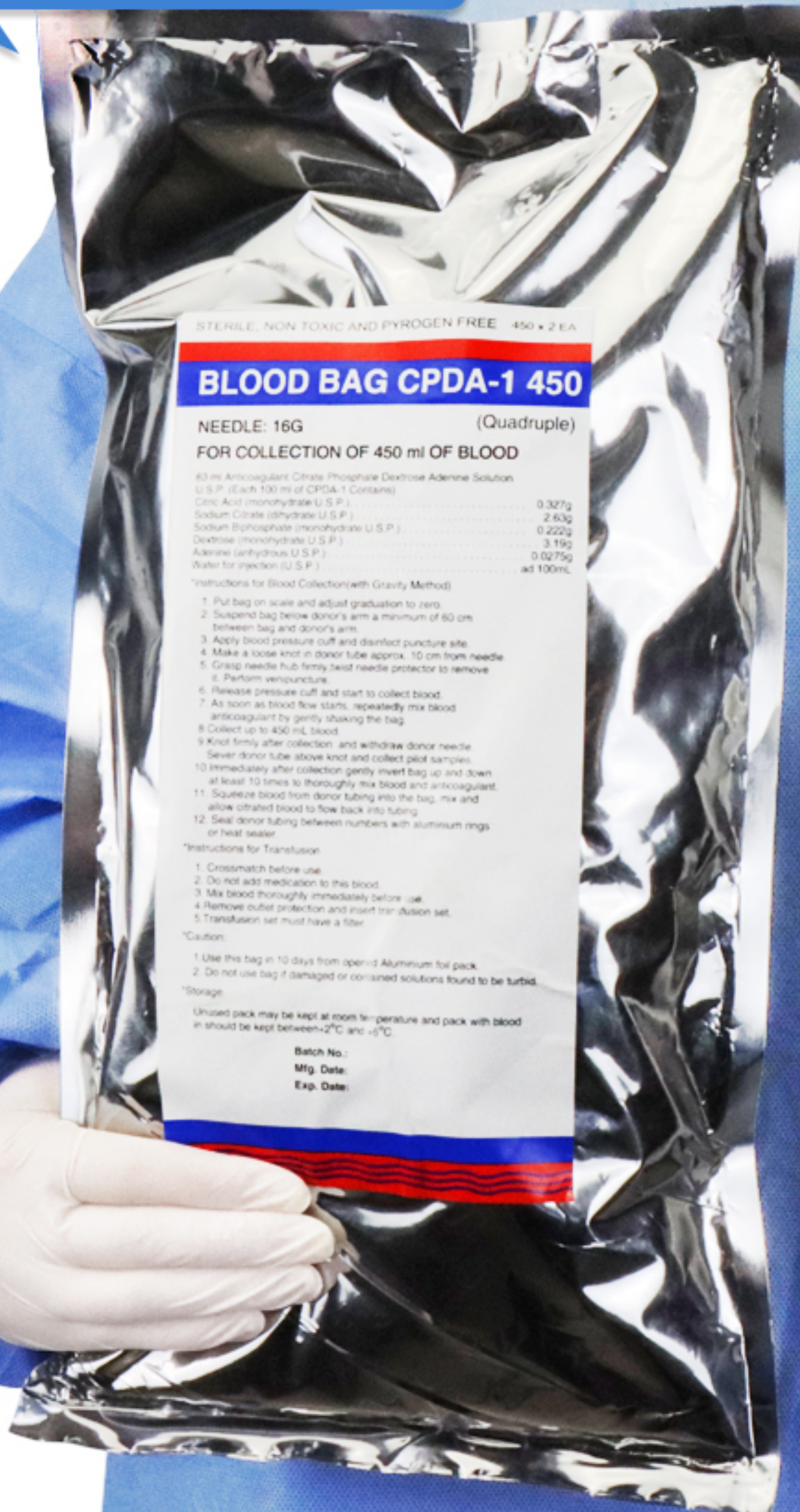

INDIVIDUAL PACKAGING



PRODUCT DESCRIPTION



BAG BODY

The body is made up of medical grade PVC materials through plastic blowing blow plastic Convenient hanger slits and holes

COLLECTING NEEDLE

16G sharp pointed Japanese needle that is siliconised with ultra-thin wall.17G needle is also available



TUBE

Provide with standard donor tubing and code number on the tube surface