

INSULIN PEN NEEDLE

Stainless steel needles imported from France



4 Sterile Needles Layer Sealed Packaging

The needles are all sterile and individually packaged, making them easy to carry around. Double protective caps on the needle seat and needle tip protect the needle body from external damage.



Sealed dialysis paper



**Needle seat
protective cover**



**Needle tip
protective cover**



Needle seat

Various Specification Available




Needle Size		Colour Identification	
34G×4mm	0.18×4mm	With	
33G×4mm	0.20×4mm	Pink	
32G×4mm	0.23×4mm	Light-green	
31G×4mm	0.25×4mm	Orange	
31G×5mm	0.25×5mm	Pink	
31G×6mm	0.25×6mm	Deep Blue	
31G×8mm	0.25×8mm	Orange	
30G×8mm	0.30×8mm	Yellow	
30G×10mm	0.30×10mm	Yellow	
29G×12mm	0.33×12mm	Red	
28G×12mm	0.36×12mm	Green	Blue
27G×12mm	0.40×12mm	Middle Grey	
26G×12mm	0.45×12mm	Brown	

Universal Connector

Suitable for all insulin injection pens.




Bevel Cut for Easy Piercing of Skin

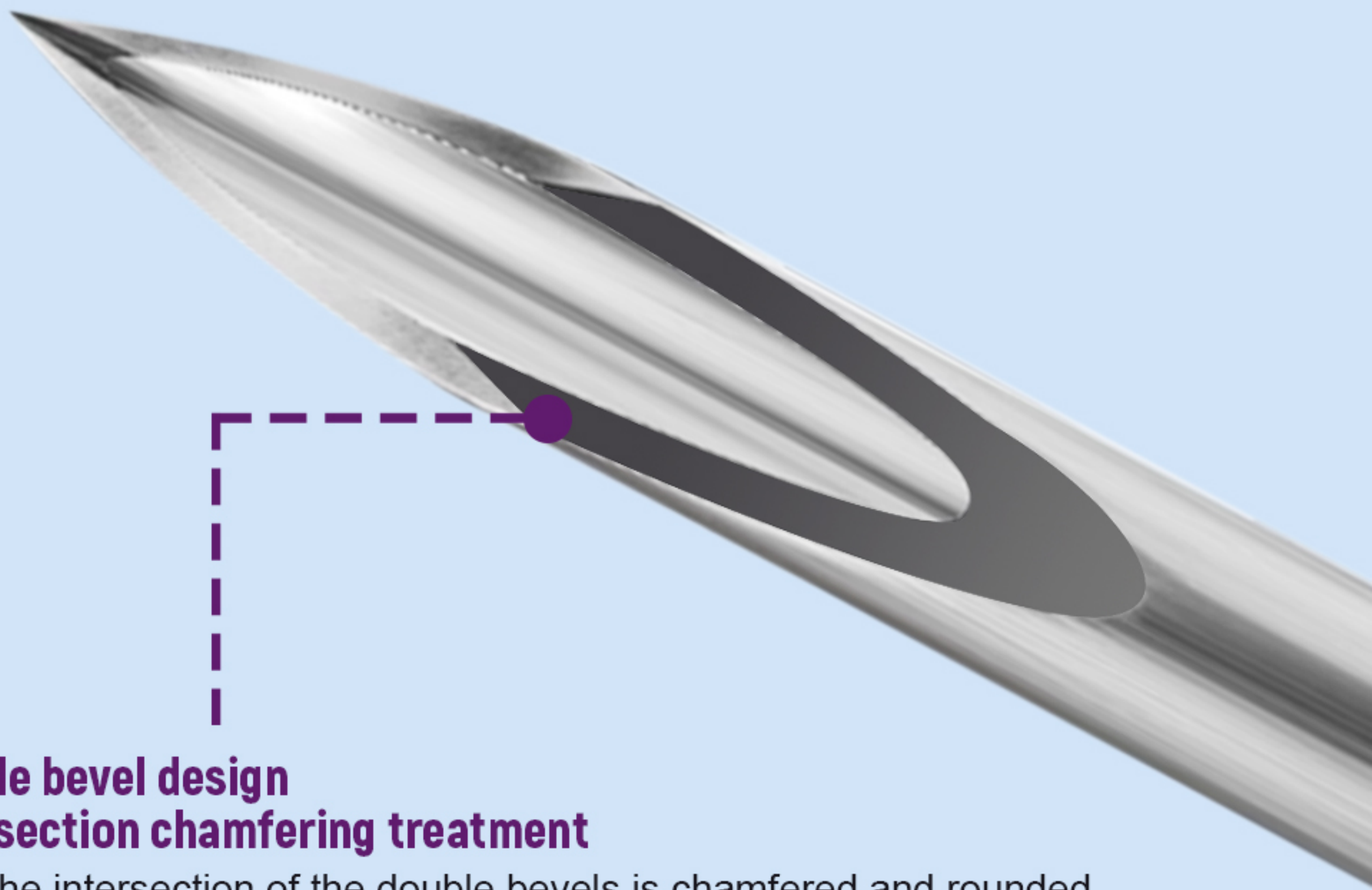

Double bevel design intersection chamfering treatment

After the intersection of the double bevels is chamfered and rounded, the needle insertion resistance is reduced and the pain of the user is reduced.

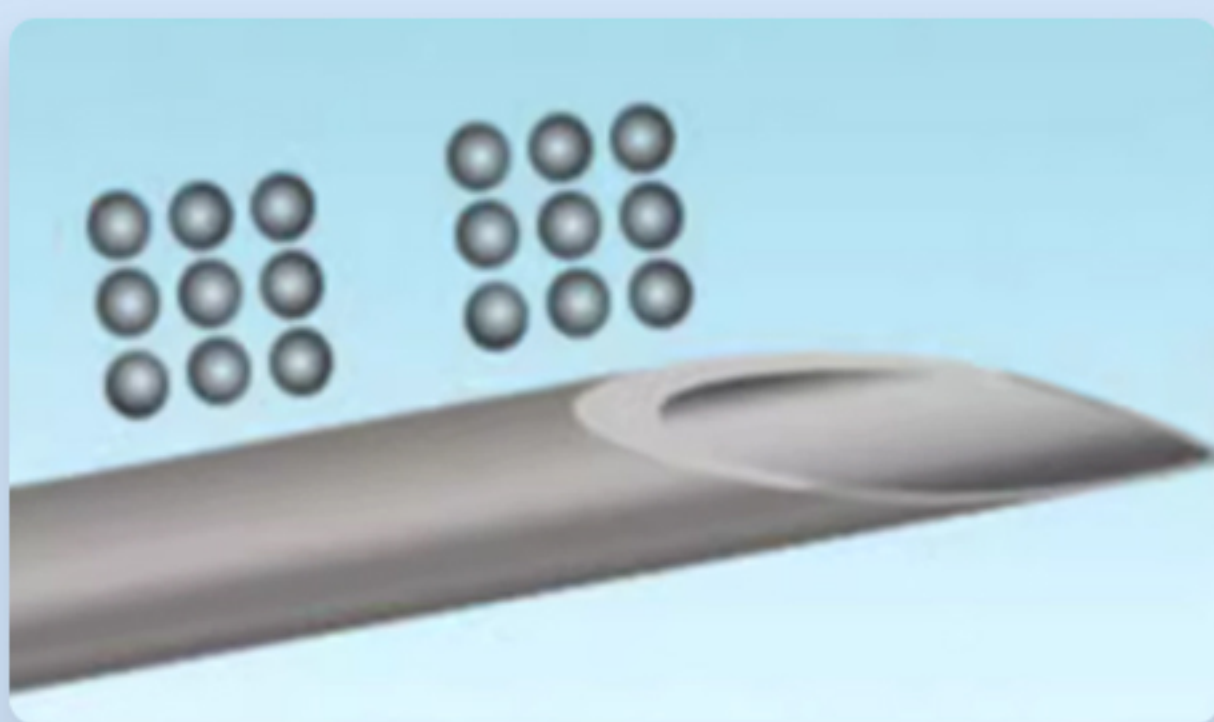

Sandblasting electronic polishing



Double siliconized lubrication treatment

